

Pipetting range 5 µL – 1.95 mL Fluid handling:

Typical reaction volume range of 200 µL to 250 µL

Probes: Liquid sensing probes monitor remaining volumes, crash protection

Accuracy of $\pm 1\% + 0.01$ A or better up to 1.5 A and ELISA reader:

±2% or better above 1.5 A

ELISA Optical filter range: 4 standard filters 405, 450, 492, 630 nm, optional additional filters

Glow Luminescence CLIA reader:

Sensitivity/Detection limit: HRP 1x10-18 moles, ALP 1x10-21 moles

Detector: Photomultiplier (PMT)

Spectral Response Range: 300-650 nm

Peak Wavelength: 400 nm

Standard sample rack: 90 sample capacity - serum vials or 13x100 mm tubes

Standard reagent rack: Typically 29 - Assorted replaceable racks and custom designed

racks are available for various bottle sizes

Integrated barcode scanning with positive patient ID system Barcode scanner:

Light source: Tungsten-Xenon lamp

Height - 57.7 cm Dimensions:

Width - 107 cm Depth - 62.4 cm







Two-Plate ELISA and Chemiluminescence

The Ultimate in Flexibility

AWARENESS TECHNOLOGY, INC.

www.awaretech.com

Automated Two-plate ELISA/CLIA Analyzer

ChemWell 2 Automated ELISA/Chemiluminescent analyzer offers a true two-plate solution for processing microwell ELISA and CLIA assays in a robust walk away system.

- Separate reagent and sample pipetting improves throughput
- Individual strips are timed separately allowing a wider variety of protocols to run simultaneously
- ChemWell 2 management software keeps track of assay timing requirements and automatically regulates assay start times to avoid conflicts. Mix and match ELISA or CLIA tests by strip row or by plate



- ChemWell 2's programmable wash uses an 8-well wash head allowing for an entire strip to be washed at once with consistent timing
- Built in barcode reader tracks sample placement
- Optional liquid sensing disposable probe tips are available for sample dispensing






- ChemWell 2's combination of hardware and software is designed to optimize ELISA/CLIA testing on one or two microwell plates.
- Follow prompts to program test kits, many flexible options. Once you have created your worklist, loaded the samples and reagents, ChemWell 2 will do the rest.
- ChemWell 2 comes with a host of reporting options, providing your laboratory with critical data in an organized printout or electronic document, reports are customizable
- Ideal for programming an assortment of microwell ELISA and CLIA kits including infectious diseases, hormones, tumor markers, auto immunity and toxins.



Holds 92 disposable tips and 90 samples

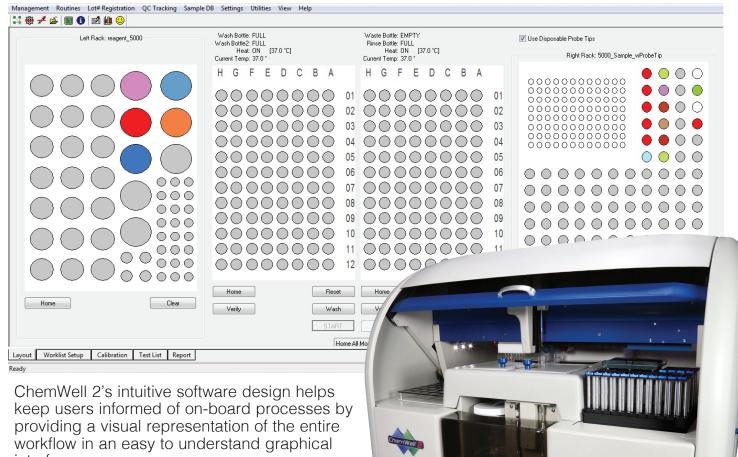

interface.

ChemWell 2's software options provide users the ultimate in *Flexibility*.